

Institute of Neuroinformatics
UNI - ETH Zurich

Arquitectura y
Tecnología de Computadores
Universidad de Sevilla

Excerpt from the Diploma Thesis
“Highspeed USB2 AER Interfaces”

USBAERmini2 and CaviarViewer Userguide

Author:
Raphael Berner

Advisors:
Dr. Tobias Delbruck (ETHZ)
Prof. Anton Civit Balcells (US)
Dr. Alejandro Linares Barranco (US)

April 14, 2006

Contents

1	Setting Up The Host Computer	1
1.1	Software setup	1
1.2	Setting the Environment Variables	1
1.3	Driver Installation	1
2	Setting Up Matlab	4
2.1	Changing The Java Classes	4
3	Setting Up The Hardware Connections	5
4	LEDs	6
5	The CaviarViewer	7
5.1	Starting and Live Displaying	7
5.2	Recording Events	8
5.3	Playing Back Recorded Files	10
6	Accessing The Java Classes From Matlab	10
6.1	Important Java Methods	11
7	Matlab Scripts	12
7.1	Important Note on Sequencing and Monitoring	13
7.2	startup.m	13
7.3	aemon.m	13
7.4	aemonseq.m	14
7.5	aeseq.m	14
7.6	aeseq_cont.m	15
7.7	aeseq_cont_stop.m	15
7.8	multi_monitor.m	15
7.9	multi_monitor_seq.m	16
8	Importing Recorded Data To Matlab	17
9	Simultaneous Sequencing with Matlab and Using CaviarViewer	17
10	Updating Firmware	17
11	Using A Device Without Firmware	18
12	Power Issues	18
13	Using USBAERmini2 Simultaneous With Other Devices	18

List of Figures

1	Select advanced driver installation.	2
2	Browse to the driver location.	2
3	Windows device manager after successful driver installation.	3
4	Netbeans Profiler Button	5
5	Netbeans Library Manager	5
6	Two USBAERmini2 boards with synchronisation connection	6
7	CaviarViewer Window, selecting an interface	7
8	CaviarViewer Window, LIVE view. Chip-type and interface are displayed in the titlebar	8
9	Reducing the driver buffer size in the CaviarViewer	9
10	Check if the time is reset in all windows.	10

List of Tables

1	USBAERmini2 operation modes	11
2	Matlab scripts	13

1 Setting Up The Host Computer

1.1 Software setup

If not yet installed, install the Java Runtime Environment 1.5. It can be downloaded from <http://www.java.com/en/download/index.jsp>.

Update or check out the code base from subversion. The URL is <https://svn.ini.unizh.ch/repos/avlasi/CAVIAR/wp5/USBAER/INI-AE-Biasgen/>.

1.2 Setting the Environment Variables

This step is not needed if you want to use the USBAERmini2 only with Matlab. But if you plan to use the CaviarViewer, it is absolutely crucial to set up the Windows PATH environment variables, so that Java can use the native code DLLs.

Open the Control Panel “System”, click on “Advanced” and then “Environment variables”. Now create a new variable called “usb2aemon”. Then open a Windows explorer window and navigate to the working copy directory of the code you checked out and copy the complete path to the two following folders, separated by a “;” (only a semicolon, no space!) to this usb2aemon variable.

```
wp5\USBAER\INI-AE-Biasgen\host\java\JNI  
wp5\USBAER\INI-AE-Biasgen\host\java\JNI\SiLabsNativeWindows
```

It will look something like this:

```
C:\Documents and Settings\tobi\My Documents\avlasi-svn\CAVIAR\wp5\  
USBAER\INI-AE-Biasgen\host\java\JNI;C:\Documents and Settings\tobi\  
My Documents\avlasi-svn\CAVIAR\wp5\USBAER\INI-AE-Biasgen\host\java\  
JNI\SiLabsNativeWindows
```

Now add this variable to the Windows PATH variable by including it surrounded by “%”, as in `...;%usb2aemon%`. PATH (or “path”, case doesn’t matter) is another environment variable. Don’t forget the semicolon “;” between the usb2aemon variable and the other PATH components but DO NOT put a space between the semicolon and the next path.

1.3 Driver Installation

After plugging in the device, the Windows New Hardware Wizard will show up. Select “No, not this time” and click on Next. Select the advanced driver installation (“Install from a list or specific location”). Select “Search for the best driver in these locations” and browse to `CAVIAR\wp5\USBAER\INI-AE-Biasgen\host\driverUSBIO`.

Windows should now be able to install the device correctly. You can check this by starting the device manager, which should now list the device as in figure 3.

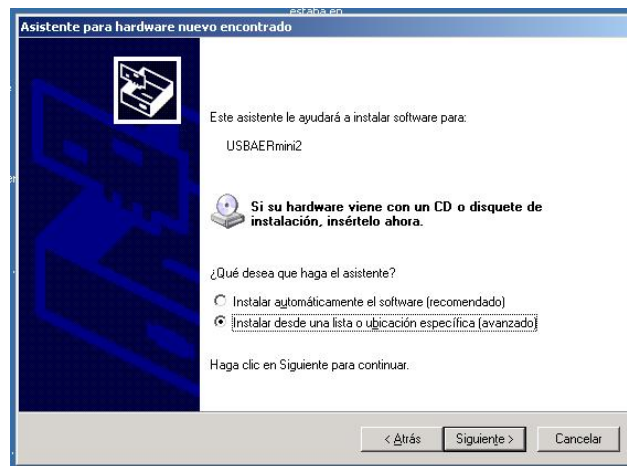


Figure 1: Select advanced driver installation.

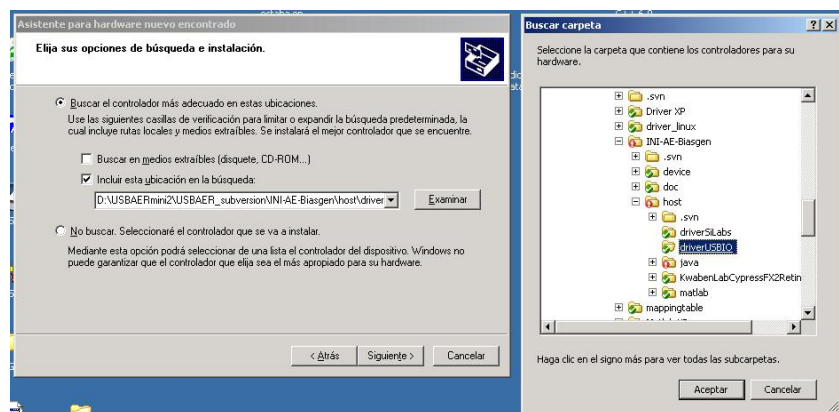


Figure 2: Browse to the driver location.

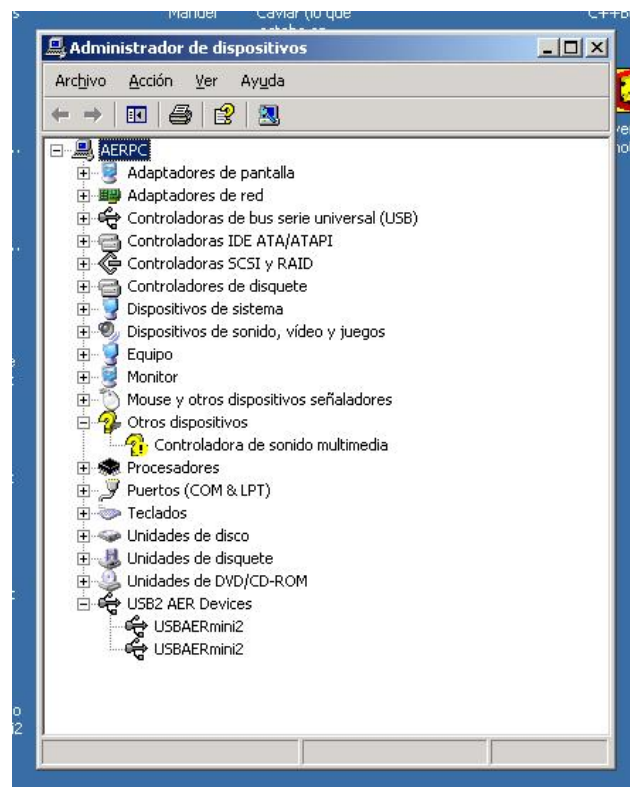


Figure 3: Windows device manager after successful driver installation.

2 Setting Up Matlab

Make sure Matlab uses Java Runtime Environment 1.5 by typing `version -java` in the command prompt. If not, create a windows environment variable called `MATLAB_JAVA` which points to the JRE 1.5 home directory.

Start matlab and edit the file `librarypath.txt` by typing `edit librarypath.txt` in the command prompt. Add the path to the folder

`CAVIAR\wp5\USBAER\INI-AE-Biasgen\host\java\JNI` to this file. If you want to use devices based on the Silicon Labs C8051F320 as well, add also the path to `CAVIAR\wp5\USBAER\INI-AE-Biasgen\host\java\JNI\SiLabsNativeWindows`. The file will then look something like this:

```
##
## FILE: librarypath.txt
##
## Entries:
##   o path_to_jnifile
##   o [alpha,glnx86,sol2,unix,win32,mac]=path_to_jnifile
##   o $matlabroot/path_to_jnifile
##   o $jre_home/path_to_jnifile
##
$matlabroot/bin/$arch
D:\USBAERmini2\USBAER_subversion\INI-AE-Biasgen\host\java\JNI
```

At this point you have to restart Matlab, otherwise the change in the `librarypath.txt` file does not become active.

Navigate to the folder `CAVIAR\wp5\USBAER\INI-AE-Biasgen\host\matlab\` and run the file `startup.m`. This adds `usb2aemon.jar` and `usbio.jar` to the dynamic matlab path and instantiates the hardware factory, through which the devices can be accessed.

Matlab is now ready to use the USBAERmini2 boards. To verify the installation, connect the sequencer output of a board directly to its monitor input and sequence a few hundred events using the script `aemonseq.m`. See section 7 for details on this script.

2.1 Changing The Java Classes

If you want to develop further functionality for the java classes, download Netbeans 5 from <http://www.netbeans.org/downloads/>. Install it and open the `usb2aemon` project directory, which is `CAVIAR\wp5\USBAER\INI-AE-Biasgen\host\java`.

You also have to download the Netbeans Profiler, which can be found at the same URL as Netbeans 5. Then you have to run the Profiler once by pressing the button shown in figure 4.

Start the Library Manager, which can be found in the *Tools* menu. You have to add two libraries, one called `usbio231`, which points to



Figure 4: Netbeans Profiler Button

INI-AE-Biasgen/host/java/jars/UsbIoJava.jar, and another called *jogl*, which points to INI-AE-Biasgen/host/java/jars/jogl.jar.

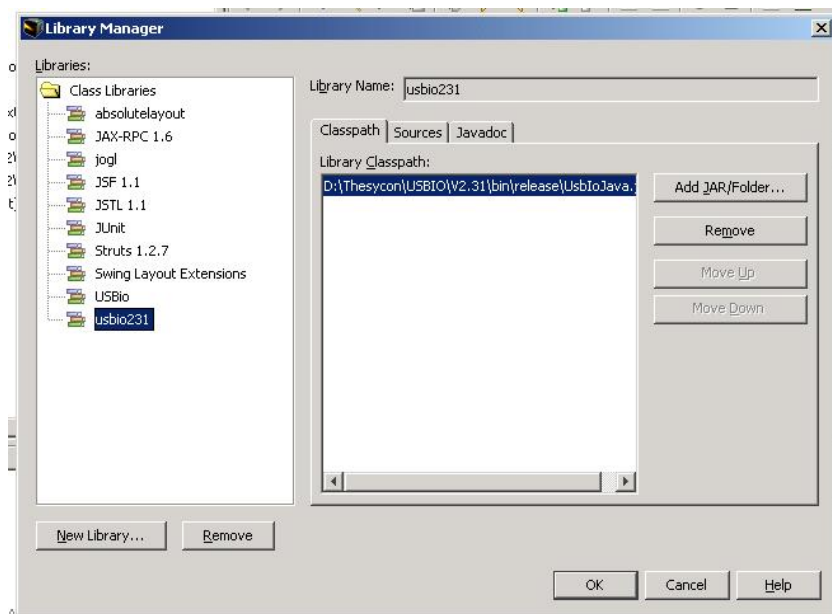


Figure 5: Netbeans Library Manager

3 Setting Up The Hardware Connections

Connect the USBAERmini2 boards to USB2.0 ports of your computer. The USBAERmini2 does not work with USB1.1. If you wish to capture synchronised data from several boards, connect the synchronisation output pin (labelled SO) from the desired master board to the synchronisation input pins (labelled SI) of the slave boards.

If the device is used in terminal mode, i.e. the monitor port is connected, but the pass-through port is not, it is best to connect request and acknowledge lines of

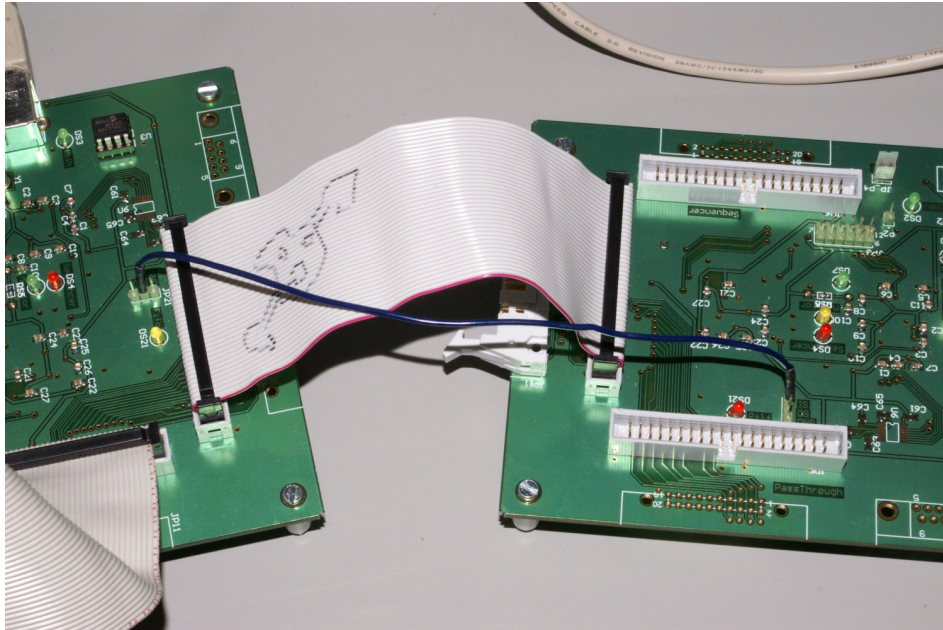


Figure 6: Two USBAERmini2 boards with synchronisation connection

the pass-through port. The device can detect if a receiver is present on the pass-through port, but only while monitoring. When monitoring is inactive, request and acknowledge lines of monitor and pass-through port are connected directly, therefore transmission is blocked if no device is present and monitoring is not active.

4 LEDs

DS1/3.3V Power LED. Should always be turned on while connected.

DS2/1.8V Power LED. Should always be turned on while connected.

DS3/BkPt Breakpoint LED. Only active in debugger mode.

DS4/Sequencer Turned on while sequencer is running in host trigger mode.

DS5/Monitor Turned on while monitor is running.

DS7 Heartbeat. Should always be blinking while connected, except while downloading to the EEPROM. If it stops, the microprocessor is stuck and the device has to be reset by unconnecting.

DS21/Fifo Master/Slave. Turned on as long as SI input is low, i.e. when the device is acting as timestamp master.

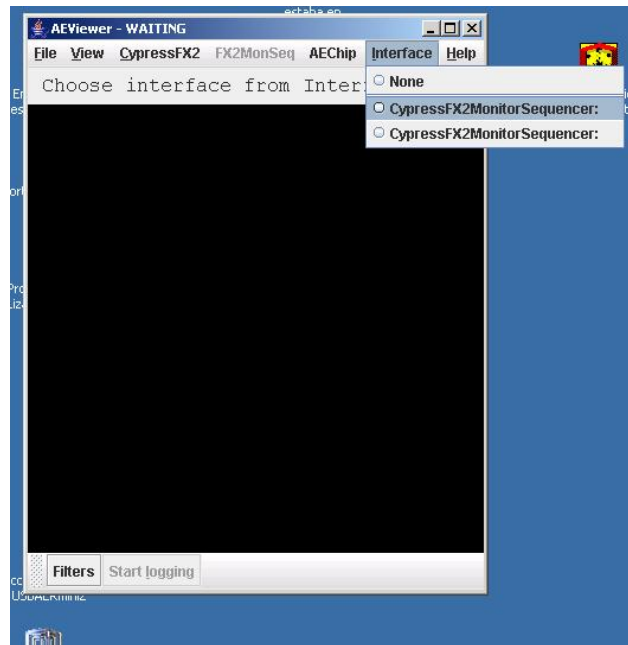


Figure 7: CaviarViewer Window, selecting an interface

5 The CaviarViewer

The CaviarViewer is a Java application to view real time data from several kinds of AER devices like the USBAERmini2, the USB2AERmapper and the TmpDiff128 retina, to record data from these devices or to view recorded data. Sequencing recorded data is not yet implemented.

5.1 Starting and Live Displaying

The CaviarViewer can be started by doubleclicking `CAVIAR\wp5\USBAER\INI-AE-Biasgen\CaviarViewer.cmd`. A window like the one shown in figure 7 should show up.

Choose the appropriate chip-type from the menu AEChip and then select the interface to display from the Interface menu, see figure 7. Note that the interfaces are sorted in the order of the PID, and secondly in the order of their name (serial number), although their name is not displayed. For example the board named *mini 1* is on top of board *mini 2*.

As soon as an interface is selected, event acquisition is started and chiptype, interface type and name are displayed in the titlebar, as can be seen in figure 8.

If the selected interface is an USBAERmini2, the FX2MonSeq menu will become active and will let you change operation mode and get the number of events missed due to full fifos. See section 6 for more details.

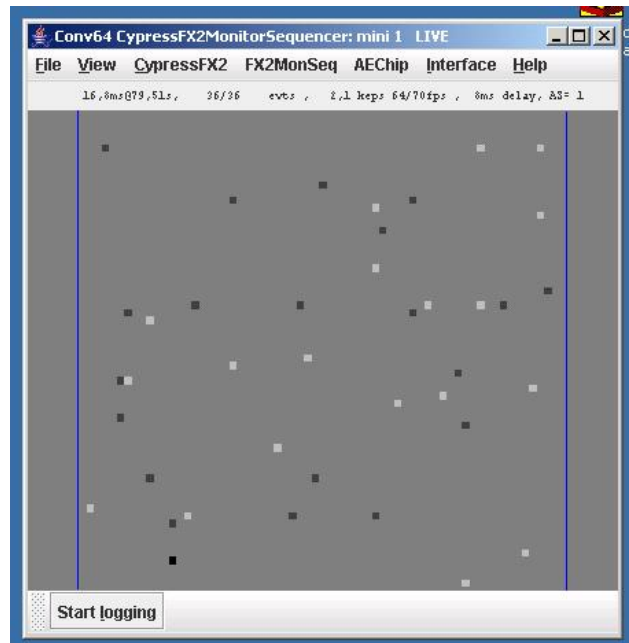


Figure 8: CaviarViewer Window, LIVE view. Chip-type and interface are displayed in the titlebar

If the display of the events seems jerky, this may be due to a low event rate, taking a long time to fill the driver’s buffers until it passes them to the application. In this case, it is best to reduce the size of the buffers. Figure 9 shows how to do this.

5.2 Recording Events

Recording of events can be started by pressing the “Start logging” button in the lower left corner and stopped by pressing it again.

Synchronized Recording of Multiple Devices

There are two types of synchronisation available.

- Software synchronisation
- Hardware (electrical) synchronisation

Software synchronisation resets the timestamps on the devices by sending a vendor request to each device in turn. This is not absolutely synchronous, but does not require the synchronisation pins to be connected. The devices are reset by pressing the zero key in one window.

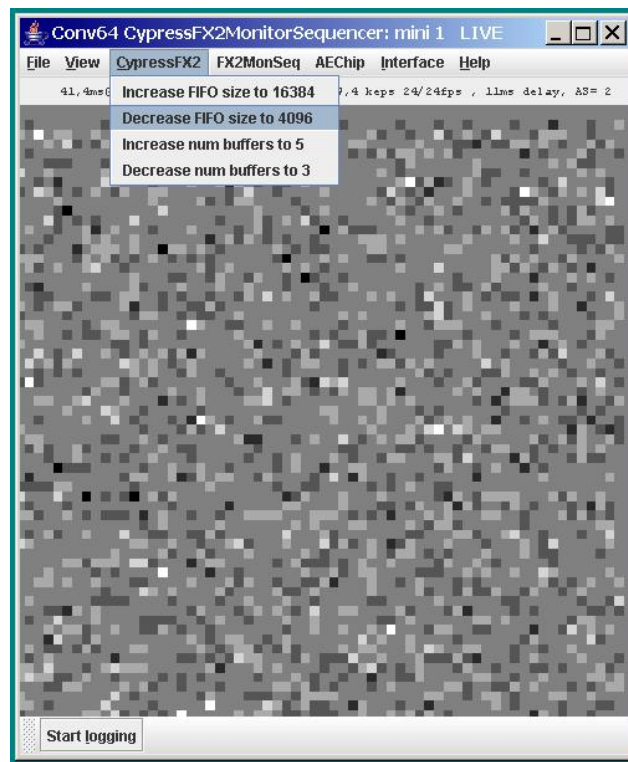


Figure 9: Reducing the driver buffer size in the CaviarViewer

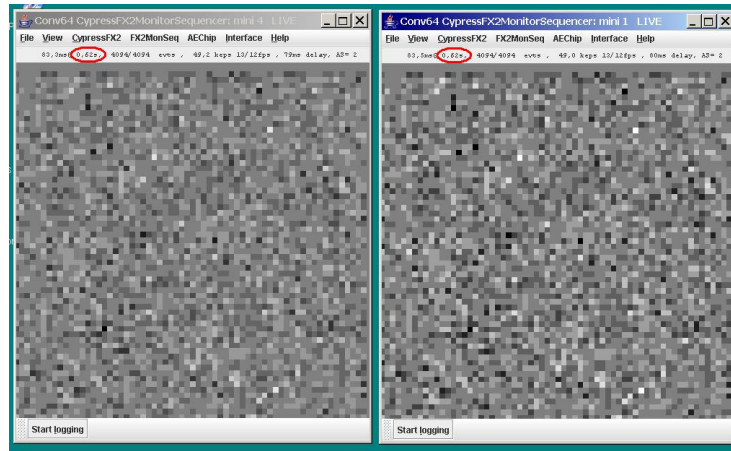


Figure 10: Check if the time is reset in all windows.

The electrical synchronisation is more accurate, because all the devices are reset from the master device. To enable electrical synchronisation, connect the SO pin of the desired master device to the SI pins of the slave devices and check the *Electrical Synchronisation* checkbox in Menu → File. Then press the zero key *in the window of the master device*. Unfortunately, the resetting does not work reliably at the moment, so check if the time really is reset in all the windows.

5.3 Playing Back Recorded Files

You can view either single data files (.dat) or synchronised sets of files (.index). You can drag and drop either type of file onto a fresh CaviarViewer AEViewer window. Or you can select the file using menu item File/Open... (shortcut "o"). If you want to select an .index file, then you need to change the file type in the file chooser; due to a bug in the graphics it only shows you one choice and doesn't show the Open button until you hover over it. This is a byproduct of using a fast native "heavyweight" Canvas to render the events.

Examine the menus for help; almost all menu items have single-key shortcuts. (e.g. r=rewind, f=faster, s=slower, etc).

6 Accessing The Java Classes From Matlab

Java classes and methods can be accessed directly from Matlab without any wrapper functions. For example the following calls instantiate a hardware factory and get a reference to the first USBAERmini2 device.

```
>> factory = ch.unizh.ini.caviar.hardwareinterface.usb.  
CypressFX2MonitorSequencerFactory.instance();  
>> usb0=factory.getInterface(0)
```

Note that the instantiation of the factory is not necessary, as this is already done by `startup.m`. See section 7 for details.

6.1 Important Java Methods

Please see the javadoc under `CAVIAR\wp5\USBAER\INI-AE-Biasgen\host\java\dist\javadoc` for more details and a complete list of methods.

CypressFX2MonitorSequencerFactory.listDevices()

List the available USBAERmini2 devices. Ignore the following errors, they are due to rapidly opening and closing the devices. *error binding to pipe for EP1 for device status: Error code 0x00000006: Windows system error code. can't set pipe parameters: Error code 0xE0001000: Operation failed.*

CypressFX2MonitorSequencerFactory.getInterface(int interfaceNumber)

This method returns a reference to a USBAERmini2 device. Parameter is the interface number, which has to be chosen from the output of `factory.listDevices`.

CypressFX2MonitorSequencer.open()

Opens the device.

CypressFX2MonitorSequencer.setOperationMode(int mode)

Sets operation mode, which includes timestamp tick and trigger mode. Valid modes are shown in table 1.

Mode	Tick	Trigger
0	1us	Master (Host)
1	0.033us	Master (Host)
2	1us	Slave
3	0.033us	Slave

Table 1: USBAERmini2 operation modes

In host trigger mode, event acquisition and sequencing is started and stopped when the commands from the host are received through USB. With this mode, it is not possible to start several devices synchronously.

In slave trigger mode, event acquisition and sequencing is started and stopped when the master device starts or stops. Therefore, slave mode works only when the SI pin is connected to the SO pin of the master device and the master device receives the start or stop command from the host. This mode is recommended to

use for the slave devices when using `multi_monitor.m` or `multi_monitor_seq.m`. These matlab scripts are described in section 7.

CypressFX2MonitorSequencer.getOperationMode()

This method returns the timestamp tick and displays the operation mode.

CypressFX2MonitorSequencer.setContinuousSequencingEnabled(bool)

Disables or enables continuous event sequencing. If enabled, the AEWriter thread rewinds at the end of the packet of events it has to send and sends these events over and over again. If disabled, the AEWriter thread send the events only once. The default state is that events are sequenced only once.

This method can be used in combination with the scripts `aemonseq.m` and `multi_monitor_seq.m`. The scripts `aeseq.m` and `aeseq_cont.m` already call this method with the suitable parameter.

CypressFX2.setAERReaderFifoSize(int size)

Sets the buffer size of the event-capturing thread. In general, a buffer size of at least 8kB leads to the highest possible event rates, however, when the CaviarViewer is used to monitor a low eventrate connection, smaller buffer sizes are needed to produce suitable frame rates.

CypressFX2MonitorSequencer.getNumMissedEvents()

This method returns an estimation of the number of events the device has missed because of the Cypress FX2 fifos being full when the host does not collect the events fast enough.

CypressFX2MonitorSequencer.writeMonitorSequencerFirmware()

This methods writes the firmware, which is saved in the `usb2aemon.jar` files as well, to the EEPROM. Use this function when a new firmware version is available.

CypressFX2MonitorSequencer.setDeviceName(String name)

Sets a new serial number. Parameter is the new serial number string. Be advised that after you plug in a device with a new serial number, the Windows New Hardware Installation Wizard will show up. The string length is limited to eight characters.

7 Matlab Scripts

Table 2 gives an overview over the matlab scripts available in `CAVIAR\wp5\USBAER\INI-AE-Biasgen\host\matlab\monitor_sequencer`.

startup.m	adds classes to the dynamic path and instantiates the hardware factory
aemon.m	monitoring a single device
aemonseq.m	monitoring and sequencing with a single device
aeseq.m	sequencing from a single device
aeseq_cont.m	continuous sequencing from a single device
aeseq_cont_stop.m	stops continuous sequencing
multi_monitor.m	monitoring from multiple devices
multi_monitor_seq.m	sequencing from one device and monitoring from multiple devices

Table 2: Matlab scripts

7.1 Important Note on Sequencing and Monitoring

Be advised that the sequencer needs *interspike intervals* for correct sequencing of spike trains, but the monitor returns *absolute* timestamps. Therefore recorded spike trains have to be processed (usually by the Matlab function `diff`) before they can be sent to the sequencer again. **IMPORTANT!**

7.2 startup.m

This script adds the path to the jar-files to the dynamic matlab path and instantiates the hardware factory, through which the devices can be accessed. This script has to be called at startup. Please note that this script is located in the parent folder, i.e. in CAVIAR\wp5\USBAER\INI-AE-Biasgen\host\matlab\.

7.3 aemon.m

```
[inaddr,ints,tick]=aemon(usbinterface,monitortime)
```

This is a script to monitor events with one USBAERmini2 device.

Parameters

usbinterface Reference to USBAERmini2 device. Get a reference to a USBAERmini2 device using for example `usb0=factory.getInterface(0)`.

monitortime How long the monitoring is active. This time is measured in seconds.

Returns

inaddr Address array returned from device.

ints Timestamps array returned from device.

tick Timestamp tick used in timestamps vector.

7.4 aemonseq.m

`[inaddr,ints,tick]=aemonseq(usbinterface,addr,ts,monitortime)`

Script to sequence and monitor events with one USBAERmini2 device.

Parameters

usbinterface Reference to USBAERmini2 device. Get a reference to a USBAERmini2 device using for example `usb0=factory.getInterface(0)`.

addr Array of addresses to be sent to device.

ts Array of interspike intervals. Note that no interspike interval should be bigger than $2^{16} - 1$. Also note that you have to set them according to the timestamp tick used on the device.

monitortime How long the monitoring is active. This time is measured in seconds.

Returns

inaddr Address array returned from device.

ints Timestamps array returned from device.

tick Timestamp tick used in timestamps vector.

7.5 aeseq.m

`aeseq(usbinterface,addresses,timestamps)`

Function to sequence events with one USBAERmini2 device. The device will stop sequencing when it has sequenced all events. Use `aeseq_cont.m` if you want to sequence continuously.

Parameters

usbinterface Reference to USBAERmini2 device. Get a reference to a USBAERmini2 device using for example `usb0=factory.getInterface(0)`.

addr Array of addresses to be sent to device.

ts Array of interspike intervals. Note that no interspike interval should be bigger than $2^{16} - 1$. Also note that you have to set them according to the timestamp tick used on the device.

7.6 aeseq_cont.m

`aeseq_cont(usbinterface,addr,ts)`

Function to continuously sequence events with one USBAERmini2 device. The device will rewind if it reaches the end of the arrays. This function is non-blocking. Call `aeseq_cont_stop.m` to stop sequencing. Use `aeseq.m` if you don't want to sequence continuously.

Parameters

usbinterface Reference to USBAERmini2 device.

addr Array of addresses to be sent to device.

ts Array of interspike intervals. Note that no interspike interval should be bigger than $2^{16} - 1$. Also note that you have to set them according to the timestamp tick used on the device.

7.7 aeseq_cont_stop.m

`aeseq_cont_stop(usbinterface)`

Stops a continuous sequencing device and releases it, so it can be accessed again from other processes.

Parameters

usbinterface Reference to the sequencing USBAERmini2 device.

7.8 multi_monitor.m

`[addr,isi,timestamps,tick]=multi_monitor(devices,monitortime)`

Function to monitor events with one or more USBAERmini2 devices.

Parameters

devices Array of references to monitor devices.

monitortime How long the monitoring is active. Time in seconds.

Returns

addr Cell array of address arrays. Same order as in monitors array.

isi Cell array of interspike intervals in nanoseconds.

timestamps Cell array of timestamps, unprocessed.

tick Timestamp tick used in timestamps arrays.

Make sure all the monitoring devices use the same timestamp tick! For synchronising the monitoring devices, connect SO pin of the desired master device to the SI pins of the slave devices, and set operation mode accordingly using `usbinterface.setOperationMode(mode)`. See 6 for more details.

7.9 multi_monitor_seq.m

```
[addr,isi,timestamps,tick]=
    multi_monitor_seq(sequencer,monitors,addr,ts,monitortime)
```

Function to sequence events with a USBAERmini2 device and monitor with one or more other USBAERmini2 devices.

Parameters

sequencer Reference to sequencer device.

monitors Array of references to monitor devices.

addr Array of addresses to be sent to device.

ts Array of interspike intervals. Note that no interspike interval should be bigger than $2^{16} - 1$. Also note that you have to set them according to the timestamp tick used on the sequencer device.

monitortime How long the monitoring is active. Time in seconds.

Returns

addr Cell array of address arrays. Same order as in monitors array.

isi Cell array of interspike intervals in nanoseconds.

timestamps Cell array of timestamps, unprocessed.

tick Timestamp tick used in timestamps arrays.

Make sure all the monitoring devices use the same timestamp tick! For synchronising the monitoring devices, connect SO pin of the desired master device to the SI pins of the slave devices, and set operation mode accordingly using `usbinterface.setOperationMode(mode)`. Use one of the monitoring devices as master and use the sequencing device in slave mode and don't forget to connect its SI pin. If you use it master mode, it will start sequencing a few milliseconds before the other devices start monitoring.

See 6 for more details on the operation modes.

8 Importing Recorded Data To Matlab

To import data into Matlab that was recorded with the CaviarViewer, use the script `loadaerdat.m` located in `CAVIAR\wp5\USBAER\INI-AE-Biasgen\host\matlab\`. A window will pop up where you can navigate to the desired file. The function returns two vectors, one containing the addresses, one containing the timestamps. Remember that the timestamp vector represents absolute timestamps, which are not suited directly for sequencing! Create a suitable vector with the following Matlab command.

```
isi = [timestamps(1); diff(timestamps)];
```

9 Simultaneous Sequencing with Matlab and Using CaviarViewer

There is one important thing that has to be aware of if simultaneous sequencing and viewing the other devices with CaviarViewer is desired.

The interface with the highest serial number should be chosen as sequencing device, because the sequencing device has to be locked while sequencing continuously, the interface menu of CaviarViewer is built only up to the interface number in front of the sequencing device. Available devices with higher serial number than the sequencer are not displayed.

Ignore the errors *CypressFX2MonitorSequencer.openUsbIo(): getDeviceDescriptor: Error code 0xE000100: Operation failed*. They are due to the driver trying to access the locked device.

10 Updating Firmware

The Cypress firmware can be updated from Matlab with the java method `CypressFX2MonitorSequencer.writeMonitorSequencerFirmware()` or from the CaviarViewer by launching the CypressFX2 EEPROM utility in the CypressFX2 menu.

To update the CPLD firmware, a Xilinx download cable and the Xilinx Impact application are required. The Impact project file is stored in the repository under `CAVIAR\wp5\USBAER\USBAERmini2\CPLD\Xilinx project\USBAERmini2\USBAERmini2.ipf`. Please note that this project file stores the location of the firmware file with an absolute path, you will have to edit it with a text editor to change the path to your checkout directory.

The current PCB unfortunately doesn't have printed information on how to connect the Xilinx JTAG download cable. The JTAG connector is JP3, the pin layout is the following:

Pin	Description
7	VCC
8	GND
9	TCK
10	TDO
11	TDI
12	TMS

11 Using A Device Without Firmware

This section describes how to handle a device without EEPROM or with an EEPROM that has not been programmed yet. This would be the case for a newly built device.

Plug the device into your computer and follow the driver installation instructions described in subsection 1.3. After successful driver installation, the device will show up as *CypressFX2Blank* in the device manager. Now start the CaviarViewer application. When the CaviarViewer application finds a blank FX2 device, it will automatically download the firmware for the TmpDiff128 retina to its RAM. Therefore you will be prompted again to install the driver. When the device is successfully installed as TmpDiff128 device, you can run the CypressFX2EEPROM utility from the CypressFX2 menu of the CaviarViewer application.

Press *Scan for device* to open the blank device, and then press *EEPROM Monitor/Sequencer firmware* to download the USBAERmini2 firmware to the device.

To download the CPLD code to the device, follow the instructions in section 10. To download the CPLD configuration, the device doesn't need FX2 firmware, but it has to be plugged into a USB port to be powered.

12 Power Issues

The USBAERmini2 was actually designed to be used with the Cypress FX2LP, which is an improved version of the FX2 and needs less power. As the FX2 draws more than 100mA from the USB bus during enumeration, some hosts may disable the device, as this is not within the USB specifications. If this is the case and the USBAERmini2 is not recognised, please try another USB port on your computer.

13 Using USBAERmini2 Simultaneous With Other Devices

Note that so far it is NOT possible to use the USBAERmini2 and the TmpDiff128 retina which uses the Cypress FX2LP simultaneous on the same computer. You must not have a retina plugged in when you want to use a USBAERmini2. Using

the USBAERmini2 with a TmpDiff128 that uses the SiLabs USB chip is no problem though.

It is also no problem to use a USBAERmini2 and a USB2AERmapper simultaneous on the same computer.